



NeoGeneStar
100 Randolph Road, Suite 2B
Somerset, NJ 08873
Tel: (732) 421-4567
Fax (908) 756-4483
www.NeoGeneStar.com

NeoGeneStar Cell Free DNA Purification kit

For purification of cell-free DNA from 10ml of plasma, serum, CSF or urine

(Additional information and protocol suggestion for cell free DNA purification is available on the website: www.NeoGeneStar.com)

For 10ml Samples, with Protease Pretreatment

For Research Use Only.

Not for human or animal therapeutic or diagnostic use.



Binding Characteristics and Sample Volume

The superparamagnetic particles bind DNA molecules from ~30 bases to > 10,000 bases.

The 10ml size NeoGeneStar™ circulating cell-free DNA Kit has been optimized for sample volumes of up to 10ml.

Catalog No	Sample Volume and Quantity	Pretreatment Buffer (20x)	NGS Protease (25X)	RNA Carrier	LYS ¹ Tubes	NGS™ Beads	Wash Buffer ²	Elution Buffer
NeoGeneStar™10ml-25-WPR	10ml x 25 preps	12.5ml	10ml	125μl	25	1.25 ml	40ml	2..5ml
NeoGeneStar™10ml-50-WPR	10ml x 50 preps	25ml	20ml	250μl	50	2.5ml	80ml	5.0ml
NeoGeneStar™10ml-100-WPR	10ml x 100 preps	50ml	40ml	500μl	100	5.0ml	160ml	10.0ml

¹LYS tubes contain chaotropic salts, which are irritants. Please wear gloves and handle with appropriate laboratory safety measures.

²Absolute ethanol must be added at 1:1 ratio prior to use for the Wash Buffer.

Procedure of the NeoGeneStar™ Circulating cfDNA Kit

For 10ml sample, add 0.5ml Pretreatment

Buffer, 0.4ml NGS Protease, and 5μl RNA

carrier to 10ml sample in a 15 ml centrifuge

tube, incubate 30 minutes at 55-60°C.

Add the entire pretreated plasma slowly

into the LYS tube, dissolve at room

temperature, then add 2ml isopropanol

and 50μl NGS™ Beads and mix well


Incubate for 30 minutes by vortexing or

inverting

Wash 2 times with Wash Buffer

Wash 2 times with 80% Ethanol

Elute

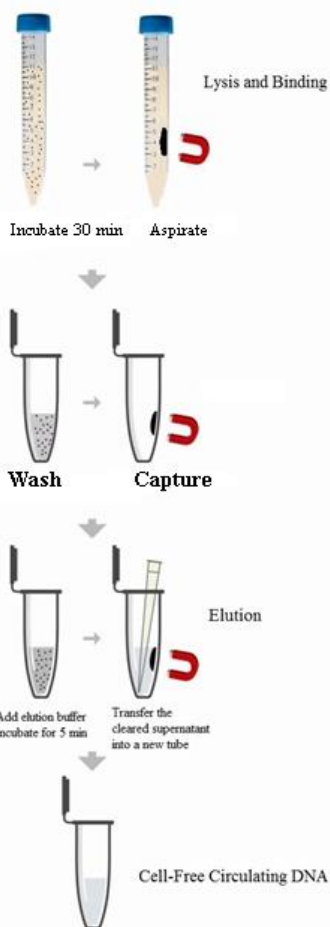



Note Regarding Wash Buffer

Wash Buffer must be diluted with an equal amount of absolute ethanol and stored at room temperature. In the event that the Wash Buffer is chilled, a precipitate may form; it would then be necessary to warm the Wash Buffer to fully dissolve the components prior to use.

Materials Needed That Are Not Supplied:

- Low DNA binding microcentrifuge tubes
- Normal Saline (0.9% w/v NaCl) solution (for volume adjustment of samples less than 10ml)
- Magnetic separation devices (for 50ml and microcentrifuge tubes)
- Absolute ethanol and isopropanol
- Tube shaker / vortexer for 50ml and microcentrifuge tubes.

Purification Protocol:

1. For 10ml sample, add 0.5ml Pretreatment buffer, 0.4ml of NGS Protease, and 5µl RNA carrier in a 50ml centrifuge tube, mix well and incubate at 55-60°C for 30 minutes.
Please note: If the sample volume is less than the designed kit processing volume, add the appropriate volume of 0.9% sodium chloride solution (normal saline) to bring the volume to the specified sample volume. (ie. if the sample volume is 9.5ml, add 0.5ml of 0.9% sodium chloride solution to result in a 10ml volume for processing sample).
2. Add the entire pretreated sample to the 10ml LYS tube and mix thoroughly at room temperature, ensure that the reagents are fully dissolved, then add 2ml isopropanol, mix well.
3. **Resuspend the NGS™ Beads by vortexing for 1 minute**, then add 50µl of NGS™ Beads, vortex or invert for 30 minutes on a vortexer or 360° rotator. Proper shaking will result in visible foam above the liquid layer.
4. Centrifuge briefly (1500rpm for 5-10 minutes) to reduce the foam, then place the 50ml tube on a magnetic stand for at least 2 minutes, the solution will clear. With the sample tube on the magnetic stand, carefully aspirate the cleared supernatant without aspirating the NGS™ Beads. A vacuum aspiration or NeoGeneStar MagPen with sleeve is convenient but not necessary.
5. Remove the sample tube from the magnetic stand, add 1ml of Wash Buffer (diluted 1:1 with absolute ethanol) into the tube and carefully rinse the NGS™ Beads and transfer to a 2.0ml low DNA binding microcentrifuge tube, then add 0.4ml diluted Wash Buffer to rinse the 50ml tube again and transfer to the microcentrifuge tube.
6. Place the microcentrifuge tube on a magnetic stand for at least 1 minute or until the solution clears, turn the tube over a few times while on the magnetic stand to remove the



NGSTM Beads trapped on the lid. With the sample tube on the magnetic stand, carefully aspirate the cleared supernatant without aspirating the NGSTM Beads.

7. Remove the microcentrifuge tube from the magnetic stand and add 1ml of Wash Buffer (diluted 1:1 with absolute ethanol) into the tube and resuspend the NGSTM Beads by vortexing for 30 seconds.
8. Place the sample tube onto the magnetic stand for at least 1 minute or until the solution clears, turn the tube over a few times while on the magnetic stand to remove the NGSTM Beads trapped on the lid. With the sample tube on the magnetic stand, carefully aspirate the cleared supernatant without aspirating the NGSTM Beads.
9. Remove the microcentrifuge tube from the magnetic stand and add 1ml of 80% ethanol to the tube. Completely resuspend the NGSTM Beads by vortexing for 30 seconds.
10. Place the microcentrifuge tube on the magnetic stand for at least 1 minute or until the solution clears, turn the tube over a few times while on the magnetic stand to remove the NGSTM Beads trapped on the lid. With the sample tube on the magnetic stand, carefully aspirate the cleared supernatant without aspirating the NGSTM Beads.
11. Wash the NGSTM Beads with 1ml of 80% ethanol again by repeating steps 9 and 10.
12. Pulse down and aspirate again to remove as much of the liquid as possible. Keeping the tube on the magnet, air-dry the NGSTM Beads at room temperature for 10 minutes.
13. Add 50µl-100µl Elution Buffer or elution buffer of your choice and resuspend the NGSTM Beads by pipetting up and down, then incubate 10 minutes at room temperature.
14. Place sample tubes on a magnetic stand, the solution will clear in about 1 minute. Transfer the cleared supernatant into a low DNA binding microcentrifuge tube. This is the purified cfDNA. Freeze the eluted cfDNA until you are ready for your downstream analysis.

Note: For certain applications, such as digital (microdroplet) PCR, it may be advantageous to briefly microcentrifuge to remove any NGSTM Beads that maybe present in the eluate.

The information in this guide is subject to change without notice.

Patent pending.

DISCLAIMER

NeoGeneStar LLC and/or affiliate(s) disclaim all warranties with respect to this document, expressed or implied. Including, but not limited to, those of merchantability, fitness for a particular purpose, or non-infringement to the extent allowed by law. In no event shall NeoGeneStar and/or its affiliate(s) be liable whether in contract, tort, warranty, or under any statute or on any other basis for special, incidental, indirect, punitive, multiple or consequential damages in connection with or arising from this document, including but not limited to the use thereof.

Notice to Purchaser: Limited use label License: Research Use Only.

The purchase of this product conveys to the purchaser the limited, non-transferable right to use the purchased amount of the product only to perform internal research for the sole benefit of the purchaser. No right to resell this product or any of its components is conveyed expressly, by implication, or by estoppel.

This product is for internal research purposes only and is not for use in commercial applications of any kind, including, without limitation, quality control and commercial services such as reporting the results of purchaser's activities for a fee or other form of consideration. For information on obtaining additional rights, please contact NeoGeneStar LLC at 100 Randolph Road, Suite 2B, Somerset, NJ, 08873. Tel: (732) 421-4567